

MultiQuant™ MD 3.0.3 HotFix 2 Release Notes



Introduction

To view information about a previous software version, refer to the *Release Notes* that came with that version of the software.

Fixed Issues

The MultiQuant™ MD 3.0.3 HotFix 2 includes:

- The fixes that were included in the MultiQuant™ MD 3.0.3 HotFix 1 for the following defects:
 - Performance issues were observed if multiple projects (>50) were stored on a network drive. (MQ-4622)
 - The addition of samples to quantitation methods generated an error. (MQ-4938)
 - The expected ion ratios were not calculated correctly after standards were removed from the Results Table. (MQ-5013)
 - Color indicators were not updated after modifications were made to the Results Table. (MQ-5321)
 - In a generated report, the value of the **Analyte Integration Type** was also shown for the **IS Integration Type**. (MQ-5597)
- Fixes for the following additional defects:
 - On workstations configured with the Windows 64-bit operating system, when an **E-Signature** was required to create or modify a quantitation method (qmethod) using the MultiQuant™ MD Software, and the electronic signature step was canceled, then the task was incomplete, and the event was not recorded in the audit trail. However, the MultiQuant™ MD Software incorrectly saved the changes to the method. (MQ-7157)
 - Under certain conditions when linear through zero regression was used, the r values were miscalculated and very small numbers or zero were shown. (MQ-7203)

Install the HotFix

| Prerequisites |
|--|
| <ul style="list-style-type: none">The MultiQuant™ MD 3.0.3 software must be installed. |

The HotFix installation folder includes installers for both the Windows 7, 32- and 64-bit operating systems and Windows 10, 64-bit operating system.

If the MultiQuant™ MD 3.0.3 HotFix 1 was previously installed, then it is automatically removed during the installation of HotFix 2. HotFix 2 includes the defect fixes from the MultiQuant™ MD 3.0.3 HotFix 1.

1. Log on to the computer as a user with Administrator privileges.
2. Make sure that the MultiQuant™ MD 3.0.3 Software is closed.
3. Make sure that the Analyst® MD Software is closed.
4. Insert the MultiQuant™ MD 3.0.3 HotFix 2 installation DVD in the DVD drive.
5. Navigate to the appropriate installation folder and then double-click the **setup.exe** file.

| Operating System | Installer |
|--|---|
| Windows 7, 32-bit system | setup.exe file located in the x86 folder |
| Windows 7 or Windows 10, 64-bit system | setup.exe file located in the x64 folder |

6. Follow the on-screen instructions to complete the installation.

Updated Files

The HotFix updates the following files in the C:\Program Files\AB SCIEX\MultiQuant 3.0\bin folder:

- Analyst32Proxy.exe (for 64-bit systems only)
- Clearcore2.ProjectUtilities.dll
- Clearcore2.RawXYProcessing.dll
- FileChecker.dll
- Sciex.MultiQuant.Data.dll
- Sciex.MultiQuant.Outliers.dll
- Sciex.MultiQuant.UI.dll

Remove the HotFix

1. Log on to the computer as a user with Administrator privileges.
2. Make sure that the MultiQuant™ MD Software is closed.
3. Make sure that the Analyst® MD Software is closed.
4. Browse to the list of installed **Programs** using the **Start > Control Panel** function and then click **MultiQuant™ MD 3.0.3 HotFix 2**.
5. Click **Uninstall**.
6. Follow the on-screen instructions.

The **MultiQuant™ MD 3.0.3 HotFix 2** is removed from the list of installed **Programs**. Removing this HotFix restores the original MultiQuant™ MD 3.0.3 Software files.

Known Issues

A slight discrepancy might be observed between data processed on computers that have different versions of .NET Framework

Data processed on computers that have .NET Framework 4.5 installed and data processed on computers that have a later version of .NET Framework installed might differ slightly. The Windows operating system automatically updates the .NET Framework to the latest version. When users reprocess data that was acquired using version 4.5 on computers with a later .NET Framework version, a slight discrepancy might be observed. (MQ 5228)

The Analyst® MD Software cannot write a checksum to the wiff file if the file remains open in the MultiQuant™ MD or other compatible software

If the Analyst® MD Software is acquiring data to a data (wiff) file, then do not open that file in the MultiQuant™ MD Software, or any other application, until acquisition has completed. If the file is opened before acquisition completes, then the Analyst® MD Software might fail to write the data file checksum. (MQ-321)

Contact Us

Customer Training

- In North America: NA.CustomerTraining@sciex.com
- In Europe: Europe.CustomerTraining@sciex.com
- Outside the EU and North America, visit sciex.com/education for contact information.

Online Learning Center

- [SCIEX University™](#)

SCIEX Support

SCIEX and its representatives maintain a staff of fully-trained service and technical specialists located throughout the world. They can answer questions about the system or any technical issues that might arise. For more information, visit the SCIEX website at sciex.com or contact us in one of the following ways:

- sciex.com/about-us/contact-us
- sciex.com/request-support

CyberSecurity

For the latest guidance on cybersecurity for SCIEX products, visit sciex.com/productsecurity.

Documentation

This version of the document supercedes all previous versions of this document.

To view this document electronically, Adobe Acrobat Reader is required. To download the latest version, go to <https://get.adobe.com/reader>.

To find software product documentation, refer to the release notes or software installation guide that comes with the software.

Note: To request a free, printed version of this document, contact sciex.com/about-us/contact-us.

This document is provided to customers who have purchased SCIEX equipment to use in the operation of such SCIEX equipment. This document is copyright protected and any reproduction of this document or any part of this document is strictly prohibited, except as SCIEX may authorize in writing.



Software that may be described in this document is furnished under a license agreement. It is against the law to copy, modify, or distribute the software on any medium, except as specifically allowed in the license agreement. Furthermore, the license agreement may prohibit the software from being disassembled, reverse engineered, or decompiled for any purpose. Warranties are as stated therein.

Portions of this document may make reference to other manufacturers and/or their products, which may contain parts whose names are registered as trademarks and/or function as trademarks of their respective owners. Any such use is intended only to designate those manufacturers' products as supplied by SCIEX for incorporation into its equipment and does not imply any right and/or license to use or permit others to use such manufacturers' and/or their product names as trademarks.



SCIEX warranties are limited to those express warranties provided at the time of sale or license of its products and are the sole and exclusive representations, warranties, and obligations of SCIEX. SCIEX makes no other warranty of any kind whatsoever, expressed or implied, including without limitation, warranties of merchantability or fitness for a particular purpose, whether arising from a statute or otherwise in law or from a course of dealing or usage of trade, all of which are expressly disclaimed, and assumes no responsibility or contingent liability, including indirect or consequential damages, for any use by the purchaser or for any adverse circumstances arising therefrom.

(GEN-IDV-09-10816-B)

For *In Vitro* Diagnostic Use.


Rx Only.

Product(s) may not be available in all countries. For more information contact your local sales representative or refer to sciex.com.

Trademarks and/or registered trademarks mentioned herein are the property of AB Sciex Pte. Ltd., or their respective owners, in the United States and/or certain other countries.

AB SCIEX™ is being used under license.

© 2020 DH Tech. Dev. Pte. Ltd.

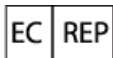

AB Sciex Netherlands B.V.
1e Tochtweg 11,
2913LN Nieuwerkerk aan den IJssel
Netherlands

MultiQuant™ MD 3.0.3 HotFix 2 Release Notes



AB Sciex Pte. Ltd.
Blk33, #04-06 Marsiling Industrial Estate Road 3
Woodlands Central Industrial Estate, Singapore 739256